



## Halogen Lamp Dental Microscopes Binocular Stereo Optical Lab Equipment A41.3406

Our Product Introduction

### Basic Information

- Place of Origin: China
- Brand Name: CNOEC, OPTO-EDU
- Certification: CE,
- Model Number: A41.3406
- Minimum Order Quantity: 1 pc
- Price: FOB \$1~1000, Depend on Order Quantity
- Packaging Details: Carton Packing, For Export Transportation
- Delivery Time: 5~20 Days
- Payment Terms: T/T, West Union, Paypal
- Supply Ability: 5000 pcs/ Month



### Product Specification

- Eyepiece Magnification: 12.5x/18B
- Adjustable Diopter:  $\pm 5D$
- Adjustable Range For Pupil Distance: 50mm~75mm
- Illumination Source: 15V/150W, Halogen Lamp
- Illumination Type: Coaxial Illumination, Illumination And Light Spot Adjustable
- Illumination:  $\geq 30000lx$
- Reaching Range Of Arm: 1230mm
- Adjustable Vertical Range: 850mm~1350mm
- Fine Focusing Speed And Range: 1mm/s~3mm/s,  $\leq 50mm$
- Horizontal Inclination Range:  $\pm 60^\circ$
- Highlight: halogen lamp dental microscopes, binocular dental microscopes, stereo dental operating microscope



### More Images



for more products please visit us on [cnoec.com](http://cnoec.com)

## Product Description

### A41.3406 Operating binocular microscope dental lab equipment

#### Product Description

#### A41.3406 OPERATION MICROSCOPE(MULTI-SECTION) dental lab equipment

---Five-step magnifications for main microscope, three-step magnifications for assistant microscope; Fine focusing for assistant microscope keeps the same clearness as main microscope.

---Lens made of imported glass from Germany can reduce the optical aberration.

---All of the lens are multi-coated, mildewproof and anti-reflection

---Equipped with two Objectives and being applicable to ophthalmology, neurosurgery, orthopaedics, microscope surgery.

--- Inclination directions and angles are controlled by hand and motor.


---Magnification and five-step fine focusing are controlled by hand and motor.

---Independent electrical box has the function of automatical diagnosis, which is very convenient for repair if necessary.

Eyepiece Magnification	12.5×/18B
Adjustable Diopter	±5D
Adjustable range for Pupil Distance	50mm~75mm
Illumination Source	15V/150W, Halogen Lamp
Illumination Type	Coaxial Illumination, illumination and light spot adjustable
Coaxial Illumination	≥30000lx
Reaching Range of Arm	1230mm
Adjustable Vertical Range	850mm~1350mm
Fine Focusing speed and range	1mm/s~3mm/s, ≤50mm
Horizontal Inclination Range	±60°
Adjustable Pitching Range(compared with vertical axis)	vertical axis) Adjusting range Depends on locking angle of the arm carrying microscope:When locked at 0°, the adjustable range -5°~+100°When locked at 30°, the adjustable range -5°~+90°When locked at 60°, the adjustable range -5°~+70°
Voltage	AC220V±10%/50Hz±1Hz, AC110V±10%/60Hz±1Hz
Power	200VA
Fuse	AC250V T3.15A, AC125V T6.3A
Electrical Safety Standard	Conform to Standard IEC 601-1. Class1 Type B
Whole Packing Volume	0.731m <sup>3</sup>
Carton Number	5
Total Weight	215kg
<b>Optional Accessory</b>	
-- Teaching tube, Beam Splitter.	
---1/3" CCD camera, CCD Adaptor, Beam Splitter and relevant Cables.	
--f400mm Objective	



Company Information

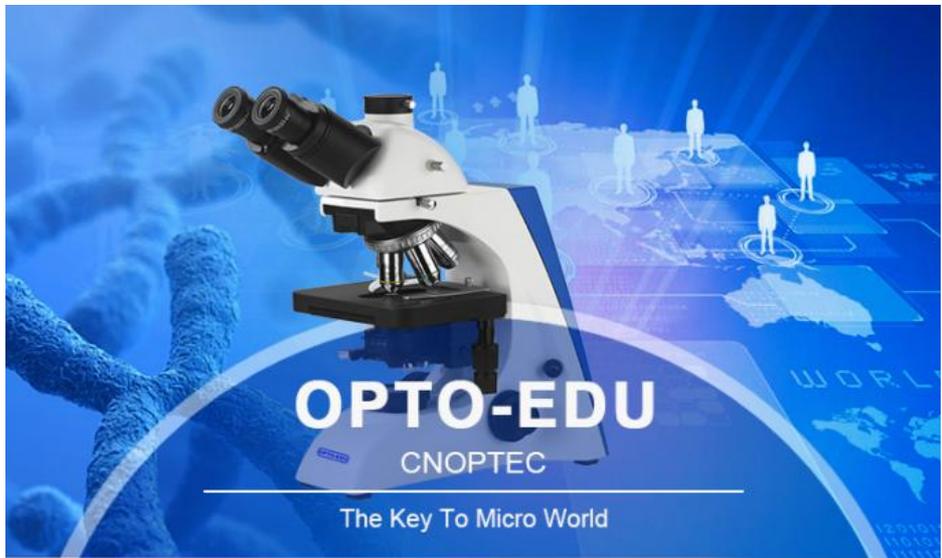

📍 OPTO-EDU is you're looking for, isn't it?

**OPTO-EDU** *cnoptec*

<p>25+ Yeas Professional Experience Know Microscopes Best</p> <p><b>150+</b> 150+ Microscope &amp; Accessory Manufacturers From</p> <p><b>25+</b></p>	<p>200+ Hot Sale &amp; Newest Update Microscope Model Every Month</p> <p><b>200+</b></p>	<p>1200+ Microscope Models To Creat One-Stop Purchase Platform</p> <p><b>720+</b> 720+ Customers From All Over World Keep Rising</p> <p><b>1200+</b></p>	<p> Top Seller For Microscope on Alibaba Recommend</p>	<p>OPTO-EDU As Famous Brand For Microscope In China &amp; USA</p> <p> Most Professional Sale Team Supply Best Microscope</p> <p></p>
---	--	--	---	---

 **Our Team And Our Factory.**

OPTO-EDU, as one of the most professional manufacturer and supplier of microscope in China, our sub-brand CNOPTEC series high end biological, laboratory, polarizing, metallurgical, fluorescene microscopes, CNCOMPARISON series forensic microscope, A63 series SEM microscope, and .49 series digital camera, LCD camera are very popular in the world market.



😊 I'm OPTO-EDU, nice to meet you!



**Opto-Edu (Beijing) Co., Ltd.**

☎ 0086 13911110627

✉ sale@optoedu.com

🌐 cnoec.com

F-1501 Wanda Plaza, No. 18 Shijingshan Road, Beijing 100043, China