



Coarse Cordless Led Light Student Compound Microscope

Our Product Introduction

Basic Information

- Place of Origin: China
- Brand Name: CNOEC, OPTO-EDU
- Certification: CE, Rohs
- Model Number: A11.1124
- Minimum Order Quantity: 1 pc
- Price: Negotiation
- Packaging Details: Carton Packing, For Export Transportation
- Delivery Time: 5~20 Days
- Payment Terms: T/T, West Union, Paypal
- Supply Ability: 5000 pcs/ Month

OPTO-EDU



Product Specification

- Eyepiece: WF10x
- Nosepiece: 3 Holes
- Objective: Achromatic 4x,10x,40x
- Working Stage: Plain Stage With Clip
- Condenser: Single Lens
- Focusing: Coarse
- Light Source: Dual Light Top & Bottom LED Light
- Power Supply: Battery
- Highlight: **led light student compound microscope , cordless student compound microscope , coarse student compound microscope**



More Images



for more products please visit us on cnoec.com

Product Description

A11.1124 Student Microscope, Cordless Dual Light, Plastic Body	
Head	Monocular Head, 360° Rotatable
Eyeiece	WF10x/16mm
Objective	Achromatic 4x,10x,40x(Spring)
Nosepiece	Triple Revolving Nosepiece
Working Stage	Plan Stage With Clips, 90x90mm
Condenser	Single Lens Condenser, With 6 Holes Disc Diaphragm
Focusing	Coarse Focusing Knob
Illumination	Dual Light Top & Bottom LED Light
Power	Cordless Powered By AA 1.5V Battery 3 pcs
Note	Order Shipped By Air Will NOT Include Batteries



**STUDENT MICROSCOPE
ECONOMICAL CORDLESS**

- 40x~400x Magnification
- Dual Up & Bottom LED Light Source
- Cordless Powered by AA 1.5V Battery 3 pcs
- Compact All in One Design for Easy Moving



A11.1124 is the best value entry level student microscope. With dual LED light, it can supply transmit or incident illumination for different object viewing. It is equipped with a plain stage with clip, and single lens condenser with disc diaphragm. It is powered by outlet 3 pcs AA batteries for cordless usage.

 MONO Monocular	 BINO Binocular	 TRINO Trinocular	 360° Head 360° Rotating	 LED LED	 ASBE Abbe Condenser	 INCIDENT Incident Light	 TRANSMIT Transmit Light	 KOHLER Kohler Illumination
 PLAN Plan Achromatic	 INFINITY Infinity Optical System	 DARK FIELD Dark Field	 POLARIZING Polarizing	 PHASE CONTRAST Phase Contrast	 FLUORESCENT Fluorescent	 ZOOM STEREO Zoom Stereo Lens	 PARALLEL Parallel Optical System	 COAXIAL Coaxial Illumination
 USB USB	 WIFI WIFI	 LCD LCD	 CORDLESS BATTERY Cordless Battery	 RECHARGEABLE Cordless Rechargeable	 WIDE POWER 100-240V Wide Power	 AUTO POWER OFF Auto Power Off	 3 YEARS 3 Years Warranty	 5 YEARS 5 Years Warranty

Product Details




01 HEAD

Monocular head inclined at 45°, rotatable 360°. With WF10x eyepiece.

02 NOSEPIECE

Triple revolving nosepiece with precise click stop. Equipped with 35mm achromatic CF objective 4X, 10X, 40Xs (oil).

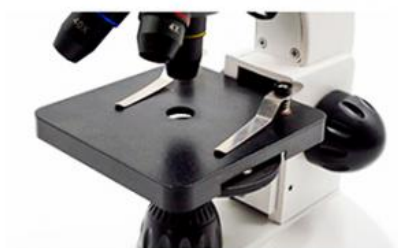

03 WORKING STAGE

Plain stage with clip, with upper safety stop to protect slide and objectives.

04 FOCUSING

Coarse focusing knob adjustment, with tension adjustment. Fine division 0.002mm, coarse stroke 36mm.



05 LIGHT SOURCE

Dual light source, with top and bottom LED light. Powered by 3 AA batteries for cordless usage.



 0086 13911110627  sale@optoedu.com  cnoec.com

F-1501 Wanda Plaza, No. 18 Shijingshan Road, Beijing 100043, China