



## 6V20W Illumination Phase Contrast Microscopy Infinity System with Plan Objective A19.0906

### Our Product Introduction

#### Basic Information

- Place of Origin: China
- Brand Name: CNOEC, OPTO-EDU
- Certification: CE,
- Model Number: A19.0906
- Minimum Order Quantity: 1 pc
- Price: FOB \$1~1000, Depend on Order Quantity
- Packaging Details: Carton Packing, For Export Transportation
- Delivery Time: 5~20 Days
- Payment Terms: T/T, West Union, Paypal
- Supply Ability: 5000 pcs/ Month



#### Product Specification

- Certification: CE
- Drawtube: Binocular
- Illumination: 6V20W Halogen Lamp
- Theory: Phase Contrast Microscope
- Working Stage: 150\*140mm
- Filter: Blue
- Highlight: **phase contrast light microscope, phase contrast microscopes**



## Product Description

### 6V20W Illumination Phase Contrast Microscopy Infinity System with Plan Objective A19.0906

- Compact Design For Appearance And Structure, with Infinity Color Correct Optical System
- Ergonomic Design, Easy Operation, Outstanding Optical Performance, Clear and Bright Image
- Specially Developed for Phase Contrast Observation, High Contrast, High Resolution, Highlight the Cell Outline and Internal Structure
- High Precision Focusing Device and Convenient Illumination Adjustment

A19.0906 Phase Contrast Microscope	
Head	Binocular, 30° Inclined, 360° Rotatable, Interpupillary Range 54~75mm
Eyepiece	PL10x/18mm, Wide Field, High Eyepoint
Phase Contrast Set	Plan Phase Contrast Objective 10X/0.25
	Plan Phase Contrast Objective 20X/0.50
	Plan Phase Contrast Objective 40X/0.65
	Plan Phase Contrast Objective 100X/1.25
	Centering Telescope Eyepiece, Dia. 23.2mm
Nosepiece	Phase Contrast Slide, for 10x, 40x Objective
	Phase Contrast Slide, for 20x,100x Objective
Nosepiece	Reversed Quadruple
Focusing	Coaxial Coarse & Fine Focusing, Low Position Design, Coarse Focusing Range 28mm, Fine Focusing Precision 0.002mm, With Tension Adjustment and Up Limit
Working Stage	Double Layer Mechanical Stage, 150*140mm, Moving Range 76*50mm, Precision 0.1mm, Double Slide Holer, Low Position Coaxial Regulating Handles
Condenser	Abbe N.A.1.25 Condenser, With Iris Diaphragm, With Filter Holder
Filter	Blue, Dia. 45mm
Light Source	Critical Illumination, 6V20W Halogen Lamp, Center Pre-adjusted
Power	90~240V Wide Range Power Supply



A19.0906 Phase Contrast Microscope






A19.0906 Phase Contrast Microscope, Side View, Front View, Back View



Binocular Head, Diopter Adjustable

#### Phase Contrast Observation

To insert a simple phase contrast flapper into the slot of condenser, then depending on the phase contrast objective, push the matching phase contrast ring into the optical path, phase contrast can be observed.

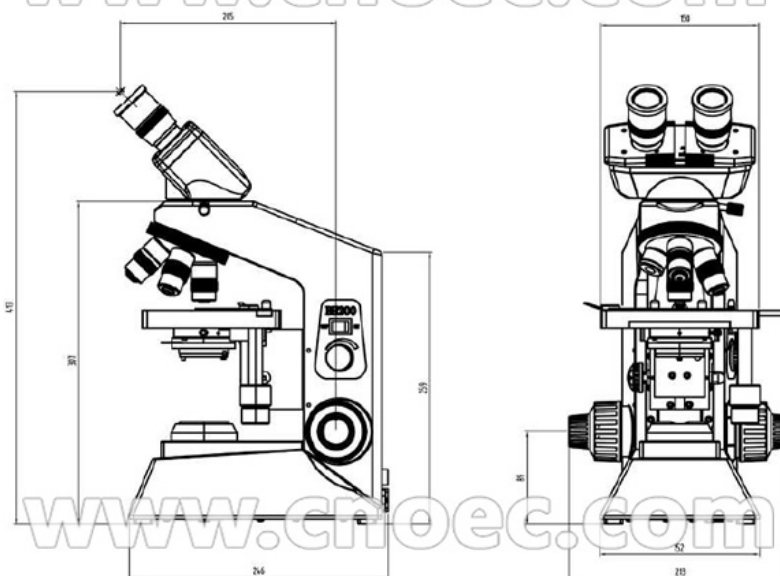

A19.0906 Phase Contrast Microscope System Size



Opto-Edu (Beijing) Co., Ltd.

☎ 0086 13911110627

✉ sale@optoedu.com

🌐 cnoec.com

F-1501 Wanda Plaza, No. 18 Shijingshan Road, Beijing 100043, China